



The Next Generation has arrived!

The all-new **Action²_{NG}** has been developed from the tried-and-tested **Action³** concept. The specifications meet the basic needs of the users whilst refusing to compromise on safety or comfort. The next generation **Action²_{NG}** shares the same essential platform as the **Action³**, offering a low total weight for easy handling combined with good driving abilities. Due to the shared use of common components, the **Action²_{NG}** perfectly fits into the **Action_{NG}** Family.



Yes, you can.®


Features and Options

Light and comfortable

With a weight from 16,2 kg, the **Invacare Action²_{NG}** fulfills all essential functional requirements. The stable folding frame means the wheelchair is easy to handle and the light weight makes transportation trouble-free. Cross-compatibility with other members of the **Action_{NG}** Family ensures that the **Action²_{NG}** can be easily adapted to meet individual user needs.

Quality of life

The tried-and-tested technology of the **Action_{NG}** Family products allow the users to go wherever they desire. Cross-compatibility with **Action³_{NG}** and **Action⁴_{NG}** make this an ideal, basic, fleet product.



Castor housing

The castor housing is mounted outside to have more space for the feet when transferring of foot propelling.



Adjustable center of gravity

The newly-designed rear wheel plate can be adjusted by 50 mm in depth and propose 3 positions in height.



Anti-tippers

Provide extra safety on slopes or uneven ground.

Technical data

Self Propulsion	380 / 405 / 430 / 455 / 480 mm	400 / 430 mm	3°	460 / 485 / 510 mm	7° / 0° - 30°	430 mm	80° / 0° - 80°	330 / 460 mm
Transit	380 / 405 / 430 / 455 / 480 mm	400 / 430 mm	3°	485 / 510 mm	7° / 0° - 30°	430 mm	80° / 0° - 80°	330 / 460 mm

Self Propulsion	200 - 280 mm (n°1) 240 mm (n°2)	Seat width + 200 mm	1030 / 1080 mm	820 / 870 mm	920 / 945 / 970 mm
Transit	200 - 280 mm (n°1) 240 mm (n°2)	Seat width + 145 mm	950 mm	725 mm	945 / 970 mm

Self Propulsion	9°	870 mm	125 kg	16,2 kg	8,9 kg
Transit	9°	870 mm	125 kg	14,9 kg	12,1 kg

All product dimensions are in millimetres to apply EN 12183 standards. Our wheelchairs comply above ISO 7176-19 standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible. For more pre-sales information about this wheelchair, please consult the www.invacare.eu.com website. It can be found under the section Medium Active Wheelchairs / Action2NG / Sales and Images in the product catalogue. For information about typical user description, intended use, adjustments and transportation, please refer to the user guide.

Invacare International Sàrl
Route de Cité-Ouest 2 - 1196 Gland
Switzerland
Tel: +41 (0)22 354 60 10
Fax: +41 (0)22 354 60 11
hqeuropa@invacare.com
www.invacare.eu.com

© 2011 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Action2NG - EU - 06/2011

Frame colour




Anthracite Black



Pearl Grey

Please notice colors may vary slightly from those displayed above.

Upholstery

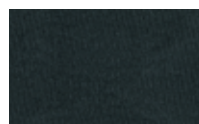

Black nylon



Yes, you can.®